



COMBI

Premium Hospital Beds



TABLES FOR WARDS

Healthcare Aids



Hospital Harmony

IMO bedside tables and overbed tables are designed with one fundamental purpose: compatibility with hospital beds, whether on a ward, ICU or other healthcare units.

With several features designed to complement the needs of bedridden patients, the combination of IMO tables and beds creates hospital environments where patients want to be and professionals like to care.

1 Combi

Epoxy painted metallic structure, with finishing options for tops and doors of lasting quality.

Overbed table with gas spring lifting option, with a lifting stroke of 960-1200 mm.

2 Adapta

Metallic structure painted in epoxy, with tops, sides and fronts in phenolic (HPL).

Overbed table with gas spring lift, with 855-1115 mm lift stroke.

3 Wally

Metallic structure, doors and drawer fronts in melamine and ABS tops.

Overbed table with gas spring lifting, with lifting stroke from 775-1025 mm.

For Any Situation

Combi support tables offer combinations for every situation. From the simple cabinet with ample storage space and a bottle rack, to the overbed table to facilitate the activities of the bedridden, or to the combined table with storage and tray.

The range as a whole is designed to meet the needs of patients and carers.



Transversal Quality

The same quality parameters applied to IMO beds are universally used in the construction of all other products and components.

This becomes apparent in the options that are available for the tables. From phenolic (HPL) to ABS plastics, all finishes have an ideal purpose and situation where they should be applied.

Cabinet

Cabinet
With Tray

Overbed
Table

Positioning

Cabinet With Tray

Rotating

Tray rotates 360° with axis in the support frame and stops at various angles.
Used by the patient to move away after use.




Lifting

Tray with 960-1200 mm lifting stroke and locking handle at desired height.
Used by the caregiver to improve patient comfort.




Extending

Tray slides up to 210 mm horizontally.
Used to compensate for bed width and centre tray with patient.

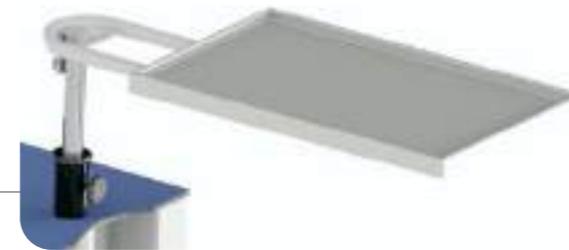

Reclining

Tray tilts bilaterally 20°, with horizontal locking.
Used by the patient to facilitate reading in bed or other activities.

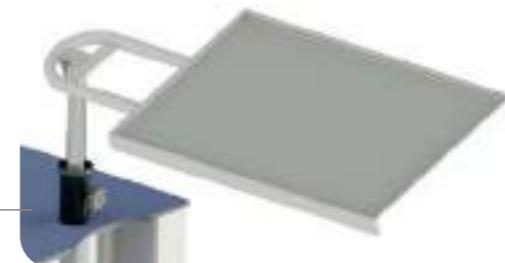

CONTINUOUS PROTECTION

Hospital Environment



Early Mobility

Simple and lightweight, the patient can easily place the cabinet against the bed by itself, as it has no handles. The smooth opening of the door and drawer to both sides helps access stored items, and the bottle holder is very easy to access.

Infection Control

Combi's simple, minimalist design makes cleaning tasks easier. While the tray and bottle holder with rounded corners prevent the accumulation of dirt, the interior features a removable shelf and plastic cuvette for thorough cleaning.

High Quality Casters

Double casters, 2 with local brake.
Option of nylon Ø50 mm or ABS Ø75 mm.



Damage Protection

Soft PVC discs on the corners to protect against accidental shocks

TECHNICAL FEATURES

Models 906.1 / 907.3

Standard - feature included as standard.
Optional - choose from various options.
Extra - feature not included as standard.

Standard	Optional	Extra
✓	✗	✗
✗	✓	✗
✗	✗	✓
✗	✗	✗

Standard - feature included as standard.
Optional - choose from various options.
Extra - feature not included as standard.

Standard	Optional	Extra
-----------------	-----------------	--------------

Structure	Robust metal frame epoxy painted	●	Base	Double ABS casters Ø75 mm, 2 local brake	●
	Doors and drawers with bilateral opening	●		Double nylon casters Ø50 mm, 2 local brake	●
	Metal doors and drawers epoxy painted	●		Corners protected by soft bumpers	●
	Phenolic (HPL) doors and drawers	●			
	Drawer with removable ABS cuvette	●	Finishes	White aluminium (doors and tops)	●
	Doors with magnetic lock	●		Signal white (structure and tray)	○
Surfaces	Table top in melamine	●		Mint Turquoise	●
	Table top in phenolic (HPL)	●		Beech	●
	Bottle rack in ABS	●		Gentiana blue	●
	Tray in melamine	●		Pigeon Blue	●
	Tray in ABS	●		Oriental Red	●
	Telescopic tray with rotation and reclining	●		Pastel Orange	●
	Manual lifting tray	●			
	Hydraulic lifting tray	●			
	Overbed table with cup tray	●			
			Registry and Conformity		

Registry and Conformity

- Certified according to the standard EN ISO 9001 - Quality Management System



CABINET WITH TRAY

* Caster Ø75mm, HPL table top

Dimensions +/- 5%

External dimensions	380 x 750 mm
Tray dimensions	375 x 595 mm
Min/Max height	960/1200* mm
Extension	210 mm
Inclination	20°



OVERBED TABLE

* Caster Ø75mm, plastic tray

Dimensions +/- 5%

Length	460 mm
Width	760 mm
Min/Max height	900/1150* mm
Inclination	20°



CABINET

* Caster Ø75mm, HPL table top

Dimensions +/- 5%

Length	565 mm
Width	380 mm
Height	730* mm


IMO - Indústrias Metalúrgicas, S.A.

Headquarters

Rua Dr. Aurélio Teixeira de Sousa, 575
4650-312 Rande Portugal
T. +351 255 340 220
sede@imo.com.pt

Showroom

Av. Defensores de Chaves, 87A
1000-115 Lisboa Portugal
T. +351 217 963 478
lisboa@imo.com.pt

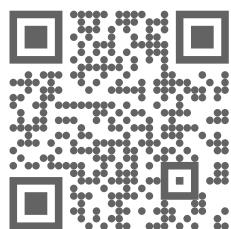

www.imo.com.pt