

FREESTYLE® COMFORT®

Portable Oxygen Concentrator with autoSAT®, UltraSense® and autoDOSE

The future is here

With its ergonomic shape, the lightweight and affordable FreeStyle® Comfort® rests against the natural curves of your body, making it comfortable to wear on the go. Equipped with user-replaceable batteries, this new oxygen therapy device meets FAA guidelines for commercial air flights and allows you to extend your time away from home.

CAIRE's proprietary autoSAT feature allows the FreeStyle Comfort to deliver a consistent volume of oxygen (bolus) per breath, up to a maximum breath rate. As the user's breath rate increases during activities like walking, the autoSAT feature adjusts the device operation to maintain a set bolus volume.



A LIGHTWEIGHT SOLUTION

Weighing only 5 lbs., the FreeStyle Comfort can go anywhere you want to go.


Enjoy the convenience of wearing it with or without your carry bag by simply attaching the no-slip, padded shoulder strap. Alternately, wear with the optional backpack.



BREATHE EASY, WORRY LESS

UltraSense, our sensitive breath detection technology, ensures oxygen is delivered to you when you need it.

Should a "no breath alarm" occur, the autoDOSE feature will pulse oxygen at 20 breaths per minute at the unit's current flow setting until the alarm state is resolved.

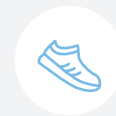

COMFORTABLE TO WEAR, EASY TO USE

The ergonomic, curved silhouette rests against the natural curves of your body, making it more comfortable to wear when you are out and about.

User-replaceable batteries allow you to extend your time away from home and stay active.



CAIRE Contact Information:

Providers Toll Free
1.800.482.2473

www.caireinc.com

SEE PRODUCT WARRANTY STATEMENT FOR COMPLETE INFORMATION. Please consult the applicable product instructions for use for product indications, contraindications, warnings, precautions, and detailed safety information.

© 2019 CAIRE Inc. All Rights Reserved.

FREESTYLE® COMFORT®

Portable Oxygen Concentrator with autoSAT®, UltraSense® and autoDOSE

Includes

- Portable Concentrator
- Carry Bag
- Battery Pack
- AC/DC Universal Power Supply and Cords
- Shoulder Strap
- UltraSense technology
- Optional autoSAT technology feature
- Other SKUs may include extra accessories

Product Features

- User-replaceable battery pack options offering up to 4 hours of battery life using the 8-cell battery pack at setting 2 and up to 8 hours of battery life with 16-cell battery pack
- Easy to wear with a shoulder strap
- Meets FAA guidelines for use on commercial air flights
- Fully functional on wall (AC) or motor vehicle (DC) power
- Easy-to-use full-color interface
- Comes equipped with the optional telehealth solution



Specifications

FLOW RATE

Pulse 1-5

DIMENSIONS

10.0 x 7.3 x 3.1 in
(25.4 x 18.5 x 7.9 cm)
8-cell battery pack

11.0 x 7.3 x 3.1 in
(27.9 x 18.5 x 7.9 cm)
16-cell battery pack

WEIGHT

5 lbs (2.3 kg)
8-cell battery pack

6 lbs (2.7 kg)
16-cell battery pack

SOUND LEVEL

40 dB(A) @ Setting 2

BATTERY DURATION

8-cell battery at setting
2 – up to 4 hours


16-cell battery at setting
2 – up to 8 hours

WARRANTY

3 years FreeStyle Comfort
(Includes sieve beds)

1 year battery and power supplies

Available Accessories

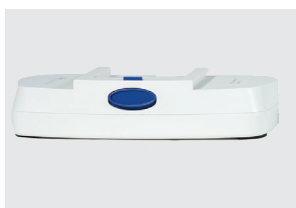

**8-CELL
BATTERY PACK**



**16-CELL
BATTERY PACK**



CARRYING BAG



BACKPACK



**DESKTOP BATTERY
PACK CHARGER
(OPTIONAL)**



**AC POWER SUPPLY
(ALL CORDS INCLUDED)**



DC POWER CORD



**EXTERNAL GROSS
PARTICLE FILTERS**

OTHER ITEMS AVAILABLE:

- Shoulder Strap
- Carry-all Accessory Bag



CAIRE Contact Information:

Providers Toll Free
1.800.482.2473

www.caireinc.com

SEE PRODUCT WARRANTY STATEMENT FOR COMPLETE INFORMATION. Please consult the applicable product instructions for use for product indications, contraindications, warnings, precautions, and detailed safety information.

© 2019 CAIRE Inc. All Rights Reserved.

ML-CONC0138 H